
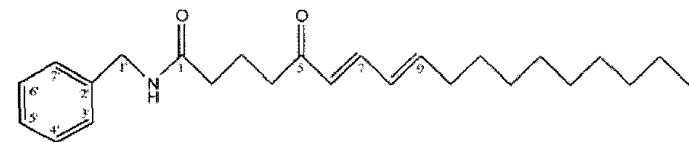



Supercritical CO2 Extracted Maca

- Purification of macamide (1)

Column: JAIGEL-ODS-AP-L, SP-120-15
Eluent : ACN
Flow rate : 5.0 ml/min
Detector : UV-310 @ 210 nm

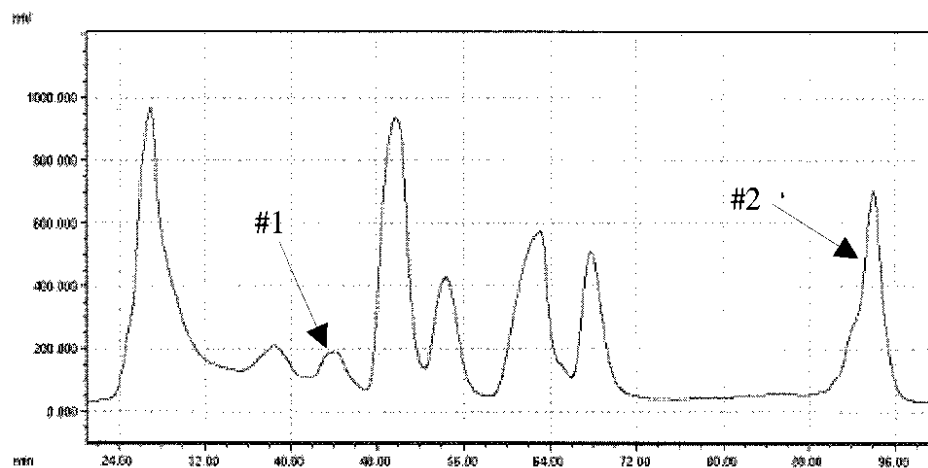

1
N-benzyl-5-oxo-6E,8E-octadecadienamide


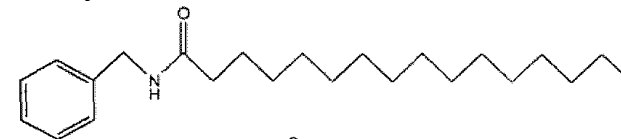


2
N-benzylhexadecanamide


Supercritical CO2 Extracted Maca

- Purification of macamide (2)

Column: JAIGEL-ODS-AP-L, SP-120-15
Eluent : ACN
Flow rate : 5.0 ml/min
Detector : UV-310 @ 210 nm

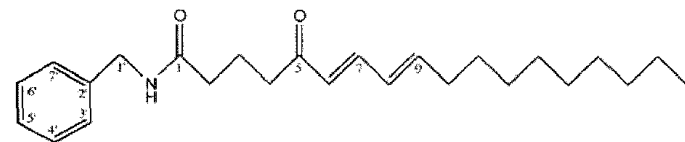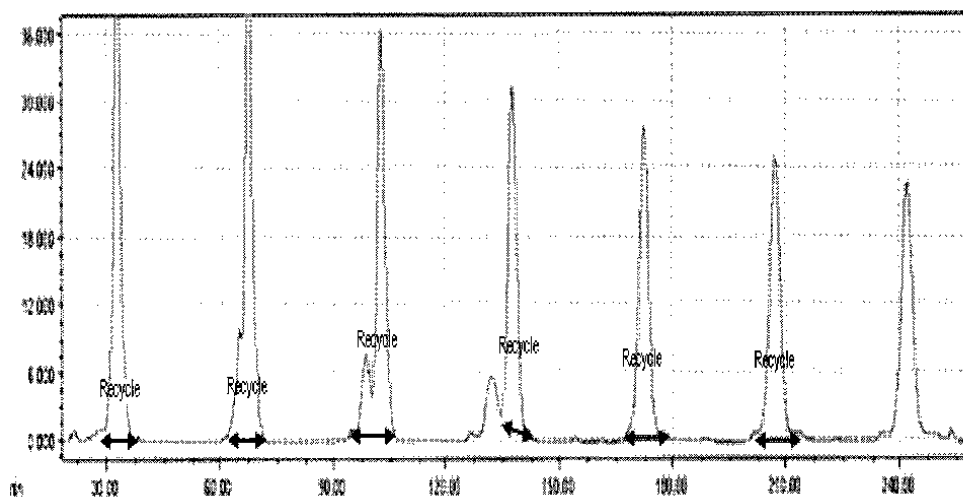

N-benzyl-5-oxo-6E,8E-octadecadienamide

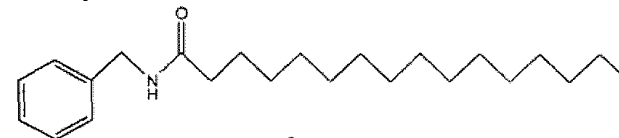

N-benzylhexadecanamide